



Lifemax

1703 - Bath Massage Pillow Instruction Manual



Lifemax Ltd, 1 Pipers Court, Thatcham, Berkshire, RG19 4ER
www.lifemaxuk.co.uk



THANK YOU FOR PURCHASING THE
**LIFEMAX BATH
MASSAGE PILLOW**

To make sure you get the maximum enjoyment from your Bath Massage Pillow, please ensure you read this instruction manual carefully.

We also recommend you keep it in a safe place should you need to refer to it in the future.

We hope you enjoy your Bath Massage Pillow.

For more products like this, please visit

www.lifemaxuk.co.uk

Lifemax **Warranty Card**
1703 - Bath Massage Pillow

In the unlikely event that your Lifemax product develops a manufacturing fault, it is covered by a one year return to base warranty when you return this card.

Your warranty is only valid with the original receipt, please ensure you keep this safe. The warranty runs from the date of purchase.

If your product develops a defect during the warranty period, please call 01635 874323 for further advice. If requested, please securely pack your item and return enclosing your contact details, your proof of purchase and fault details.

We will, at our discretion, repair or exchange the item in line with warranty legislation. Where no proof of purchase is provided, or the product is outside the warranty period, repairs/replacements will be offered at cost to the customer.

Name	
Address	
Telephone	
Email	
Retailer	
Date of Purchase	

We would like to keep you up to date with other helpful products. If you would rather we did not, please tick here.

Contents

Introduction

What does it do?

What's inside?

Functions & Controls

Before use...

How to use...

Precautions

Troubleshooting

Technical Info

Introduction

Bring even more relaxation to your bath time with our unique, IPX7 rated, waterproof massager. Featuring a carefully choreographed program, our bath massage pillow provides a gentle vibration massage – ideal for soothing painful backs, aching shoulders and tired calves.

Powered by 4 x AA batteries (not supplied), the massage pillow offers two levels of intensity delivered by three massage motors. The comfortable cover can be removed for washing and features a useful drying loop.

Compact and lightweight, the bath massage pillow is totally portable, so you can take your little piece of relaxation wherever you go.


What does it do?

- Helps you to relax your whole body whilst bathing.
- Stimulates acupuncture points to reduce fatigue and promotes physical and mental relaxation.
- Massage has been proven to improve blood circulation.
- Helps you to relax tense nerves, in turn reducing stress.
 - Helps to reduce muscle soreness and stiffness

What's inside?




Instruction Manual

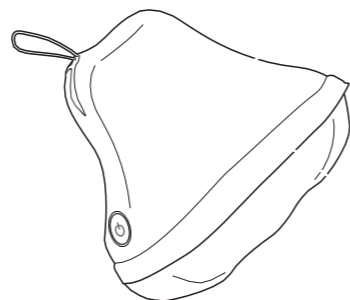


Bath Massage Pillow

Functions & Controls

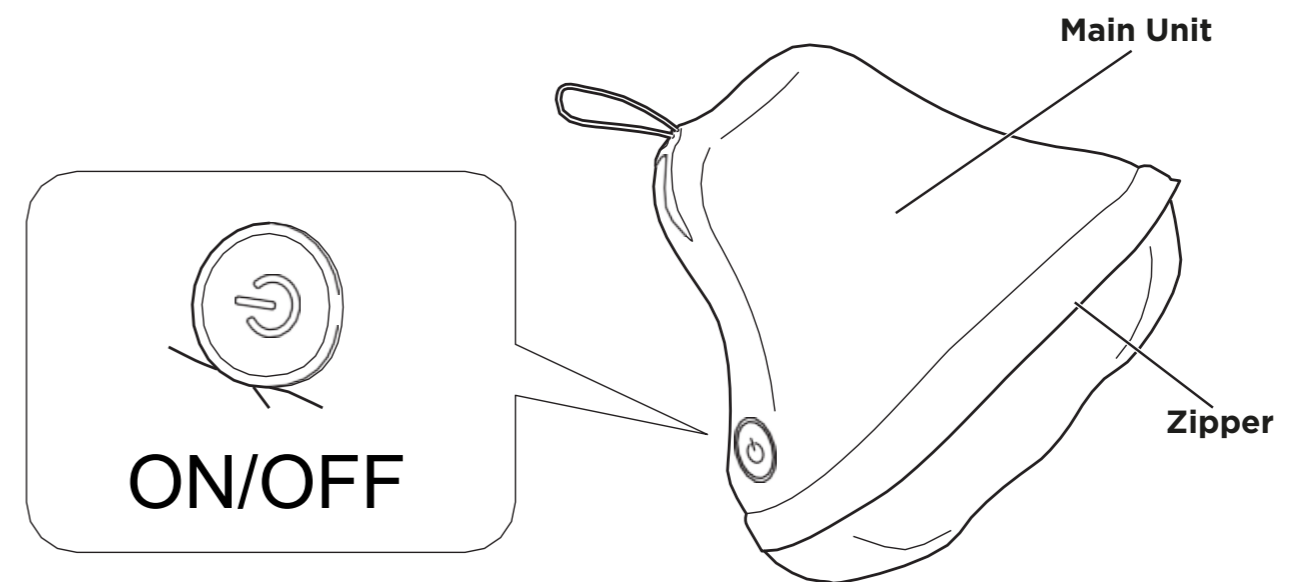

To turn on the massager, press and hold the power button.

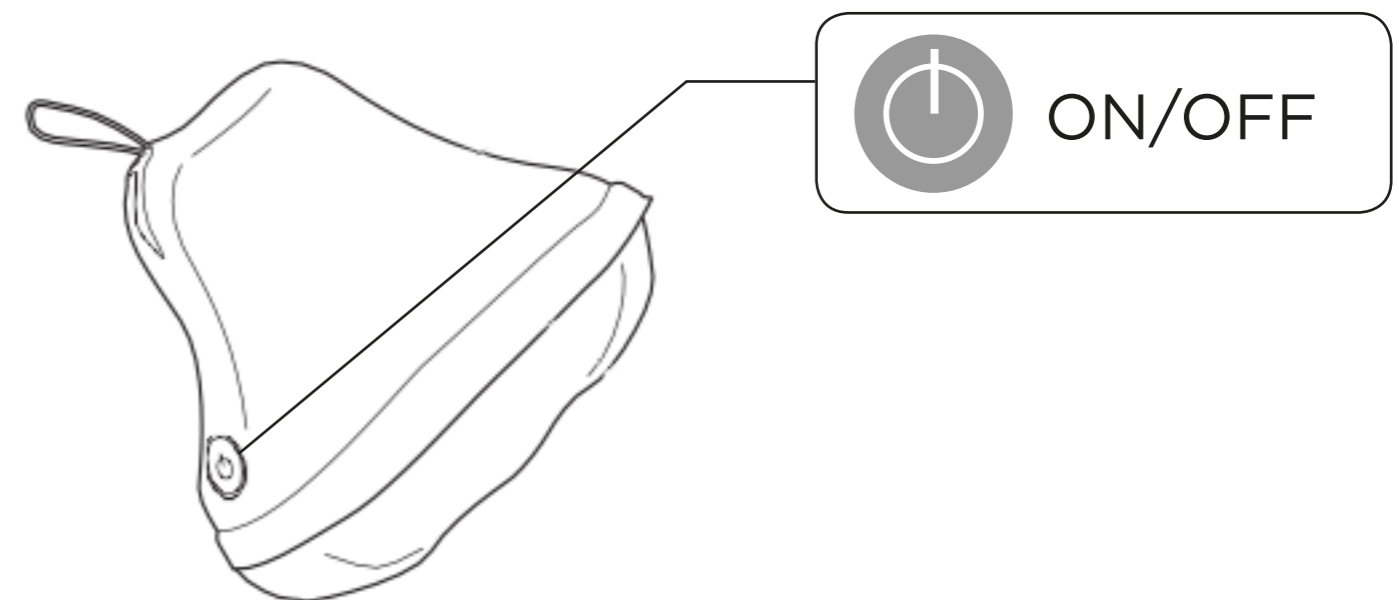


When you turn the massager on it will default to the low power massage mode. Press the power button once to switch to high power mode.


To turn the massager off, press and hold the power button.

Before use...


Install 4 x AA batteries as follows...

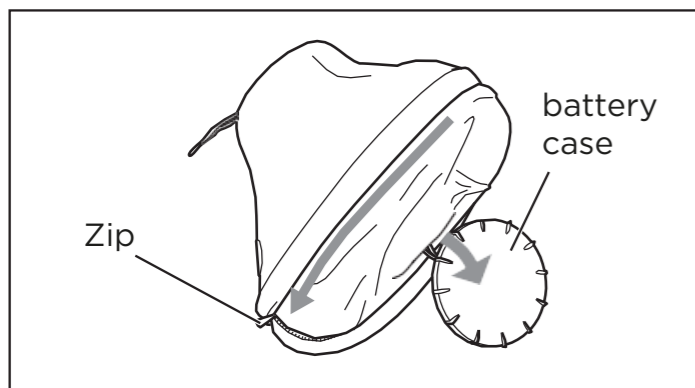


1. Take the battery case out of the massager

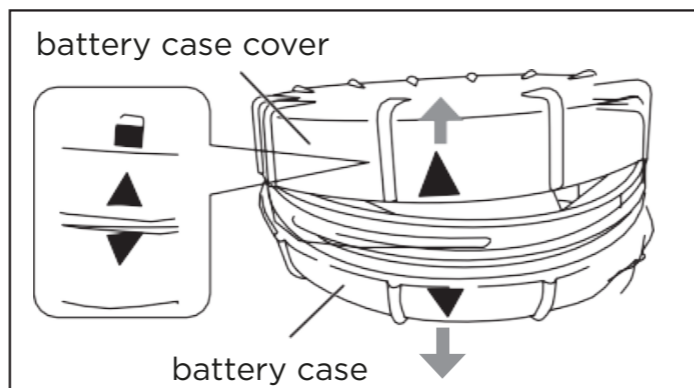


2. Open the battery case as shown above

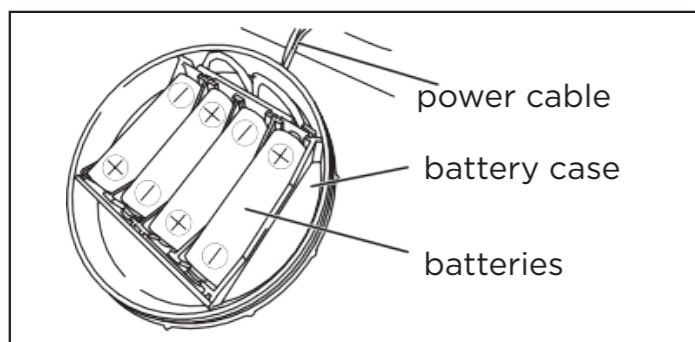


3. Insert 4 x AA batteries observing the correct polarity

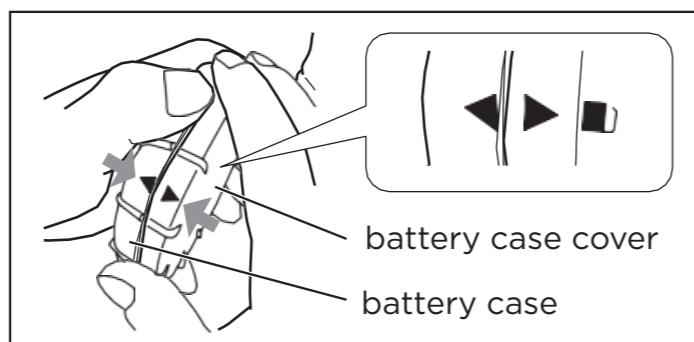

4. Close the battery case as shown above

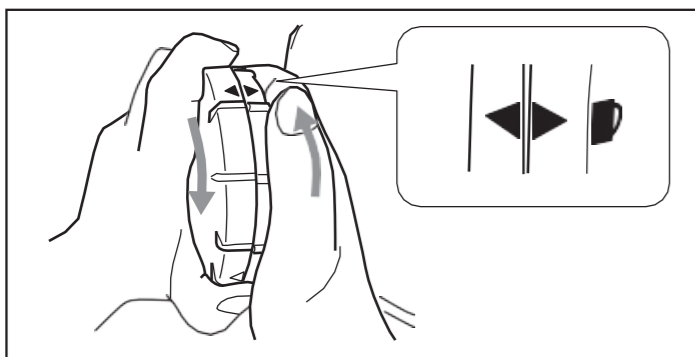

5. Check the case is closed correctly and place back into the massager

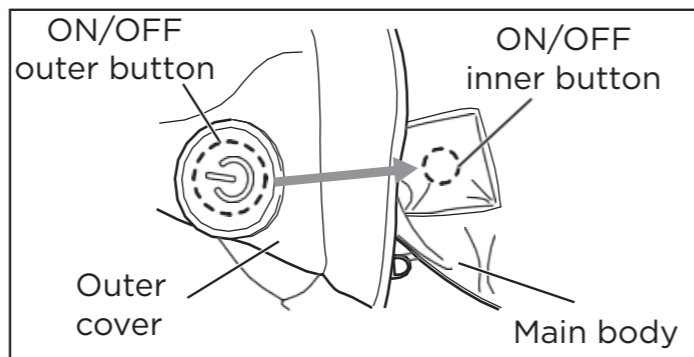

6. Ensure the inner and outer power button are aligned and zip up the cover.

How to use...

Run yourself a nice hot bath, turn your massager on and use it anywhere on your body that feels good! Simple!

Precautions

Please read this Instruction Manual carefully to ensure you are using the massager correctly. Please also note:

- Do not use whilst driving or walking
- Store the massager away from humidity, heat and/or corrosive gas.
- Individuals with symptoms such as inflammation, blood stasis or those who have undergone a recent operation should not use this product.
- Do not use foot ornaments or jewellery when using this product
- Do not apply nail polish within an hour before using this product.
- If the massager will not be used for a long time, remove the batteries to avoid leakage.
- This massager is designed as a personal care product. It is not for diagnosis / treatment of any health disorder or injury, Nor it is meant for commercial use.
- Do not machine wash the outside cover. Only wash by hand.
- Use only the correct size batteries ensuring the correct polarity is observed.

Troubleshooting

Problem	Possible Reason	Possible Solution
Massager will not power on	Batteries are low / run out	Change batteries
	Batteries are not installed correctly	Check if batteries have been inserted correctly and are seated suitably.
Weak vibration	Batteries are low / run out	Change the batteries. Check if batteries have been inserted correctly and are seated suitably.

Technical

Product Name: Bath Massage Pillow

Product SKU: 1703

Power Supply: 4 x AA batteries

Ideal Environment: Temperature: 5°C - 40°C Relative humidity: ≤80%

Material: Main body: Sponge Foam, PVC waterproof cover **Outside cover:** Refer to the washing label.

Product size: 300mm X 260mm X 80mm

Waterproof Grade: IPX7

